

AQUASOFT* Nitrile Exam Gloves



A new level of comfort, protection and value

AQUASOFT* COMFORT

AQUASOFT* Gloves are comprised of a unique nitrile formulation that provides users with both comfort and tactile sensitivity. Like our entire nitrile glove family, AQUASOFT* Nitrile Exam Gloves are easy to don with a silky, snug fit.

Key features include:

- Preferred tactile sensitivity
- Textured fingertips
- Beaded cuffs
- Ambidextrous

AQUASOFT* Nitrile Exam Gloves are also HACCP-certified and suitable for food handling, which allows facilities to simplify glove product mix ordering.¹

AQUASOFT* PROTECTION

HALYARD* AQUASOFT* Nitrile Exam Gloves provide you with the confidence of increased protection by providing resistance to the following chemotherapy drugs:

The following chemotherapy drugs and concentrations had NO breakthrough detected up to 240 minutes:

Bendamustine HCl (5 mg/ml)	Fulvestrant (50mg/ml)
Bleomycin Sulfate (15 mg/ml)	Gemcitabine (38 mg/ml)
Busulfan (6 mg/ml)	Idarubicin (1 mg/ml)
Carboplatin (10 mg/ml)	Ifosfamide (50 mg/ml)
Carfilzomib (2mg/ml)	Irinotecan (20 mg/ml)
Cetuximab (2mg/ml)	Mechlorethamine HCl (1 mg/ml)
Chloroquine (50 mg/ml)	Melphalan (5mg/ml)
Cisplatin (1 mg/ml)	Methotrexate (25 mg/ml)
Cladribine (1mg/ml)	Mitomycin C (0.5 mg/ml)
Cyclophosphamide (20 mg/ml)	Mitoxantrone (2 mg/ml)
Cyclosporin A (100 mg/ml)	Oxaliplatin (2 mg/ml)
Cytarabine (100 mg/ml)	Paclitaxel (6 mg/ml)
Cytovene (Ganciclovir) (10mg/ml)	Pemetrexed (25mg/ml)
Dacarbazine (10 mg/ml)	Retrovir (10 mg/ml)
Daunorubicin (5 mg/ml)	Rituximab (10 mg/ml)
Decitabine (5mg/ml)	Topotecan HCl (1 mg/ml)
Docetaxel (10 mg/ml)	Trisenox (Arsenic Trioxide) (1 mg/ml)
Doxorubicin HCl (2 mg/ml)	Velcade (Bortezomib) (1 mg/ml)
Epirubicin (2 mg/ml)	Vincristine Sulfate (1 mg/ml)
Etoposide (20 mg/ml)	Vidaza (Azacytidine) (25mg/ml)
Fludarabine Phosphate (25 mg/ml)	Vinorelbine Tartrate (10mg/ml)
Fluorouracil (50mg/ml)	Zoledronic Acid (0.8/ml)

The following hazardous drugs (opioids) and concentration had NO breakthrough detected up to 240 minutes:

Fentanyl Citrate injection (100 mcg/2 ml)

The following chemotherapy drugs and concentration showed breakthrough detected in less than 40 minutes:

Carmustine (3.3 mg/ml): DO NOT USE

ThioTEPA (10 mg/ml): DO NOT USE


WARNING: Not for use with Carmustine, ThioTEPA

The AQUASOFT* Powder-Free Nitrile Exam Glove material was tested with the above chemotherapy drug concentrations in accordance with ASTM D 6978, Standard Practice for Assessment Of resistance of Medical Gloves to Permeation by Chemotherapy Drugs.²



AQUASOFT* VALUE

HALYARD* AQUASOFT* Nitrile Exam Gloves are available in a 300-count dispenser box that contains three times the glove volume of a 100-count box in approximately the same dispenser box footprint.³

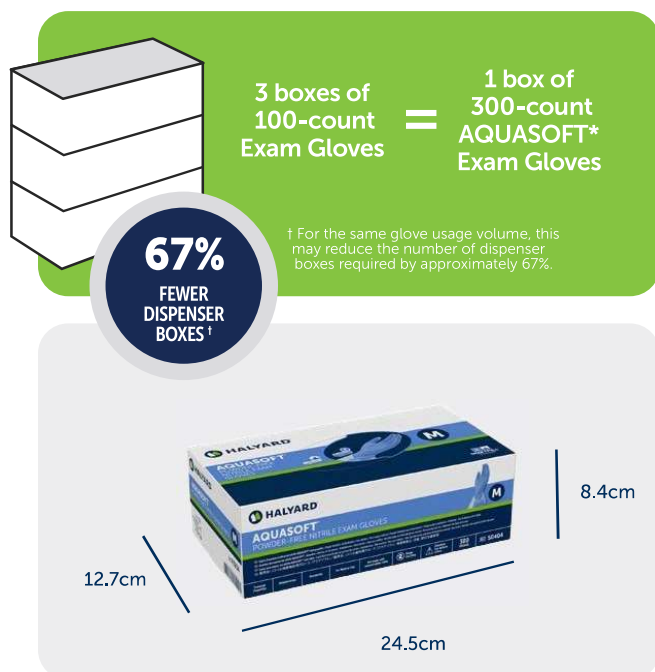

	HYH CODE	PGC/CH2 CODE	DESCRIPTION	PACKAGING
	50201	2833364	Exam Glove Dispenser, (Suits 9.5" - 12" gloves), Clear Plastic, Single, 96mm (D) x 144mm (W) x 253mm (H)	Eaches = 1
	50202	2833372	Exam Glove Dispenser, (Suits 9.5" - 12" gloves), Clear Plastic, Double, 96mm (D) x 274mm (W) x 256mm (H)	Eaches = 1
	50203	2833381	Exam Glove Dispenser, (Suits 9.5" - 12" gloves), Clear Plastic, Triple, 96mm (D) x 406mm (W) x 256mm (H)	Eaches = 1

HALYARD* AQUASOFT* Nitrile Exam Gloves

HYH CODE	PGC/CH2 ORDERING CODE	DESCRIPTION	SIZE	PACKAGING
50404	2540799	AQUASOFT* Nitrile Exam Gloves	X-SMALL	300 eaches/box 10 boxes/case
50405	2540801	AQUASOFT* Nitrile Exam Gloves	SMALL	300 eaches/box 10 boxes/case
50406	2540810	AQUASOFT* Nitrile Exam Gloves	MEDIUM	300 eaches/box 10 boxes/case
50407	2540828	AQUASOFT* Nitrile Exam Gloves	LARGE	300 eaches/box 10 boxes/case
50408	2540836	AQUASOFT* Nitrile Exam Gloves	X-LARGE	250 eaches/box 10 boxes/case

PHYSICAL PROPERTIES⁴

Before Aging

Tensile Strength:	23MPa
Ultimate Elongation:	528%

After Aging

Tensile Strength:	23MPa
Ultimate Elongation:	423%

TARGET SPECIFICATIONS⁴

Gauge Thickness Measurements	MM	MIL
Middle Finger:	0.06	2.6
Palm:	0.05	2.2
Cuff:	0.04	1.8
Average Length:	240	9.45"

QUALITY STANDARDS⁴

Exceeds current ASTM Standards for critical defects, tensile strength and puncture resistance.
Passes primary skin irritation testing in accordance with ISO 10993-10(2010).
Manufactured in accordance with Quality System ISO 13485.

AQL – ASTM D 5151⁴

Critical Defects Level	1.0 AQL [†]
[†] Acceptable Quality Level	

ASTM D6319 Standard Specification for Nitrile Examination Gloves for Medical Application. ISO-10993-10 PRIMARY Skin Irritation & Sensitization. ASTM F739 Standard Test Method for Permeation of Liquids and Gases through Protective Clothing. Materials under Conditions of Continuous Contact. ASTM F1671 Standard Test Method for Resistance of Materials Used in Protective Clothing to Penetration by Blood-Borne Pathogens Using Phi-X174 Bacteriophage Penetration as a Test System.

CAUTION: The testing conditions used are intended to approximate the worst case conditions for clinical use. Testing was conducted on single layer glove material. It is the users' responsibility to determine the applicability of these gloves for their intended use with chemotherapy drugs.

Storage Recommendations: Store in a cool, dry place. Open box should be shielded from exposure to direct sunlight, intense artificial light, x-ray machines and other sources of ozone.

O&M Halyard Australia Pty Ltd makes no representation that its products described in this catalogue will provide complete protection from exposure to contaminants, fluids or micro-organisms or will eliminate the risk of contracting infectious diseases.

**Our
Commitment
to you**

**Better
Infection
Protection**

www.halyardhealth.com.au

O&M Halyard Australia Pty Ltd
ParagonCare Customer Care Number: 1300 724 274

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE.

1. CLAIM-08250
2. CLAIM-07590
3. CLAIM-08076
4. Data on file

*Registered Trademark or Trademark of O&M Halyard or its affiliates. ©2026 O&M Halyard International. All Rights Reserved.

COPY-05901